

# **CytoFluor™**

**STAT1 (Phospho-Ser727)  
Fluorometric Cell-Based ELISA Kit**

**Catalog #: FLUO-CBP1678  
Reactivity: H:S727, M:S727, R:S727**

**Storage: 4°C  
Expiration: 6 months from receipt**

**For Research Purposes Only. Not  
Intended for Diagnostic or Clinical  
Procedures.**

**Manual Version: 23.58.4587**

# CONTENTS

# P A G E

---

|                                                          |    |
|----------------------------------------------------------|----|
| Introduction.....                                        | 3  |
| <i>Fluorometric Cell-Based ELISAs</i>                    |    |
| <i>STAT1 (Phospho-Ser727) Cell-Based ELISA</i>           |    |
| Assay Format.....                                        | 4  |
| Materials Included/ Storage and Stability.....           | 5  |
| Antibody Specificity.....                                | 6  |
| <i>Anti-STAT1 (Phospho-Ser727) Antibody Western Blot</i> |    |
| <i>Anti-STAT1 Antibody Western Blot</i>                  |    |
| <i>Anti-GAPDH Antibody Western Blot</i>                  |    |
| <i>Anti-STAT1 Antibody ELISA with Stimulants</i>         |    |
| Reagent Preparation.....                                 | 8  |
| Assay Protocol.....                                      | 9  |
| Data Analysis.....                                       | 11 |
| Short Protocol.....                                      | 12 |
| Health and Safety Precautions/Assay Restrictions.....    | 13 |
| Troubleshooting Guide.....                               | 14 |
| ELISA Plate Template.....                                | 15 |
| Technical Support.....                                   | 16 |

## INTRODUCTION

### ***Fluorometric Cell-Based ELISA***


The Fluorometric Cell-Based ELISA Kit allows for the detection of various target proteins and the effects that certain stimulation conditions have on target protein expression in different cell lines. Qualitative determination of target protein concentration is achieved by an indirect ELISA format. In essence, the target protein is captured by target-specific primary (1°) antibodies while Dye 1-conjugated and Dye 2-conjugated secondary (2°) antibodies bind the Fc region of the 1° antibody. Through this binding, the dye conjugated to the 2° antibody can emit light at a certain wavelength given proper excitation, hence allowing for a fluorometric detection method.

A monoclonal antibody specific for human GAPDH is included to serve as an internal positive control in normalizing the target RFU values. If a phosphorylated target is being detected, an antibody against the non-phosphorylated counterpart will be provided for normalization purposes. The RFU values obtained for the non-phosphorylated target can be used to normalize the RFU value for the phosphorylated target.

### ***STAT1 (Phospho-Ser727) Fluorometric Cell-Based ELISA***

The STAT1 (Phospho-Ser727) Cell-Based ELISA Kit is a convenient, lysate-free, high throughput and sensitive assay kit that can monitor STAT1 protein phosphorylation and expression profile in cells. The kit can be used for measuring the relative amounts of phosphorylated STAT1 in cultured cells as well as screening for the effects that various treatments, inhibitors (ie. siRNA or chemicals), or activators have on STAT1 phosphorylation.

## ASSAY FORMAT

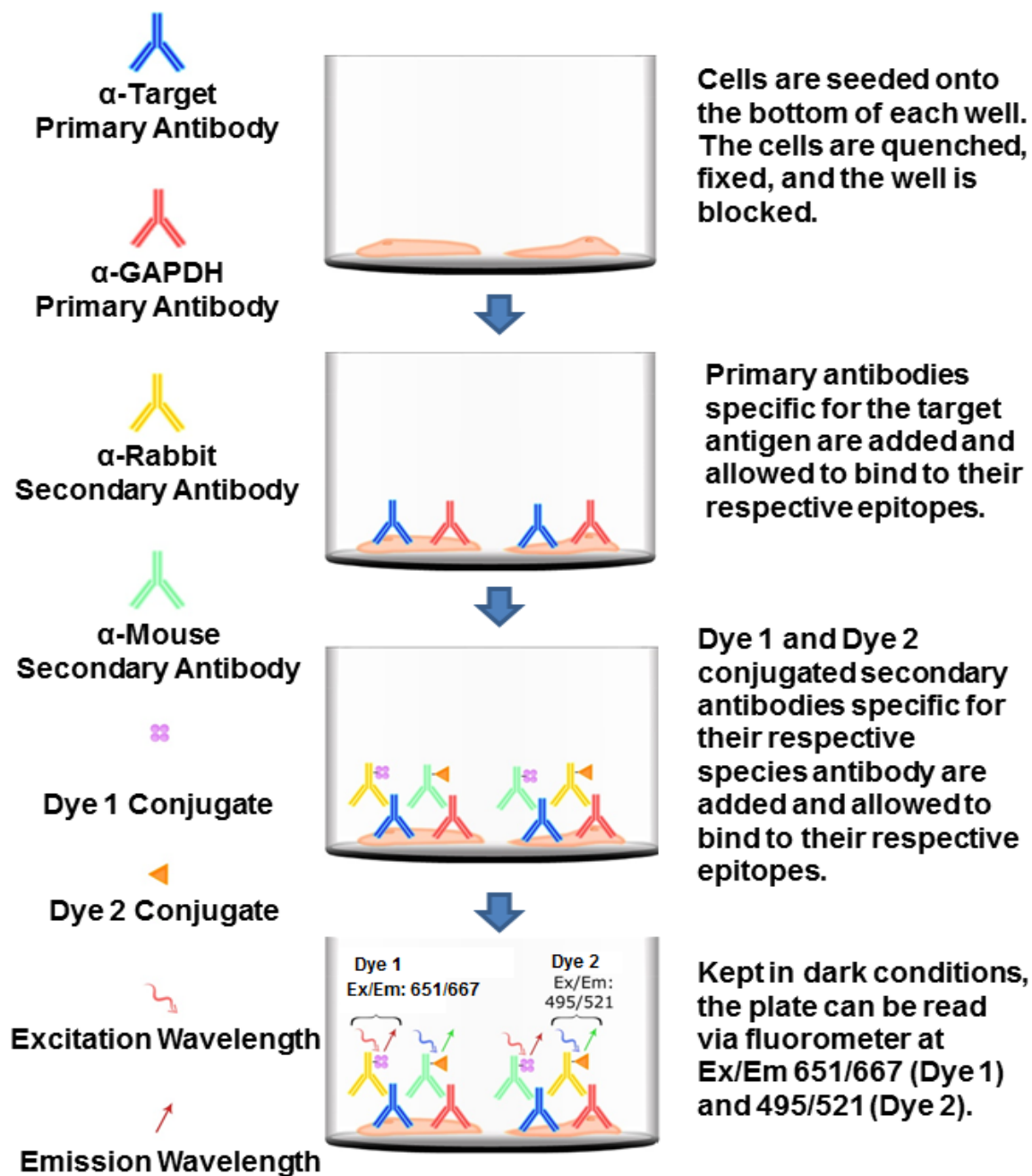

## MATERIALS INCLUDED


| Component                                          | Quantity | Storage/Stability after first use |
|----------------------------------------------------|----------|-----------------------------------|
| 10x TBS                                            | 30 ml    | 1 month at 4°C                    |
| Quenching Buffer                                   | 24 ml    |                                   |
| Blocking Buffer                                    | 50 ml    |                                   |
| 15x Wash Buffer                                    | 60 ml    |                                   |
| Primary Antibody Diluent                           | 12 ml    |                                   |
| 100x Anti-STAT1 (Phospho-Ser727) Antibody          | 60 µl    |                                   |
| 100x Anti-STAT1 Antibody                           | 60 µl    |                                   |
| 100x Anti-GAPDH Antibody                           | 110 µl   |                                   |
| Dye-1 Conjugated Anti-Rabbit IgG Antibody          | 6 ml     |                                   |
| Dye-2 Conjugated Anti-Mouse IgG Antibody           | 6 ml     |                                   |
| 96-Well Black Cell Culture Clear-Bottom Microplate | 2 Plates | -                                 |
| Adhesive Plate Seals                               | 2 Seals  | -                                 |

## ADDITIONAL MATERIALS REQUIRED

The following materials and equipment are NOT provided in this kit but are necessary to successfully conduct the experiment:


- Fluorescent plate reader with two channels at Ex/Em: 651/667 and 495/521
- Micropipettes with capability of measuring volumes ranging from 1 µl to 1 ml
- **37% paraformaldehyde stock concentration or 8% paraformaldehyde from other sources**
- Deionized or sterile water
- Squirt bottle, manifold dispenser, multichannel pipette reservoir or automated microplate washer
- Graph paper or computer software capable of generating or displaying logarithmic functions
- Absorbent papers or vacuum aspirator
- Test tubes or microfuge tubes capable of storing ≥1 ml
- Orbital shaker
- **Poly-L-Lysine or Poly-D-Lysine for suspension or loosely attaching cells**

ANTIBODY SPECIFICITY

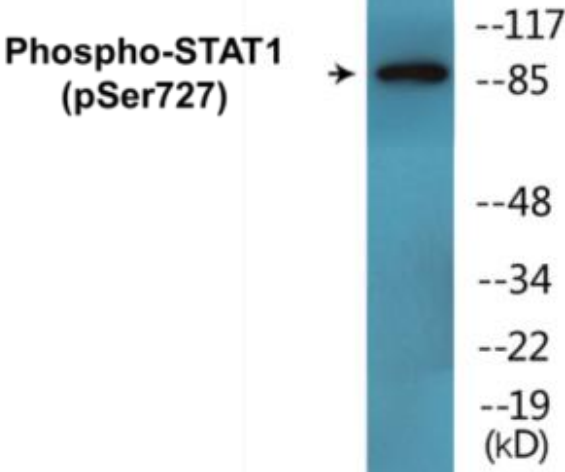


**Figure 1. Anti-STAT1 (Phospho-Ser727)**

Western blot analysis of extracts from 293 cells, 3T3 cells treated with UV (15mins) and Jurkat cells treated with eto (25uM, 24hours), using STAT1 (Phospho-Ser727) Antibody.

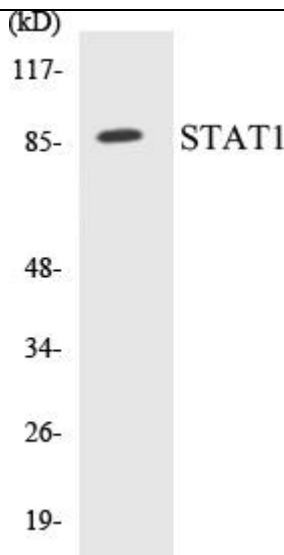

**Figure 2. Anti-STAT1 Antibody WB**

Western blot analysis of extracts from 293/3T3/Jurkat

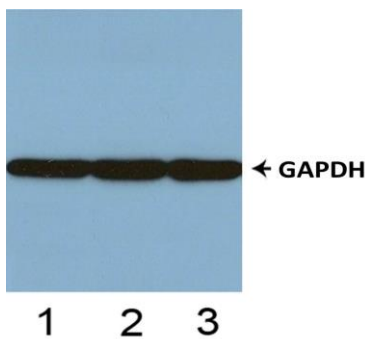


**Figure 3. Anti-GAPDH Antibody WB**

Western blot analysis of tissue lysates from human (1), mouse (2) and rat (3).

## ELISA SPECIFICITY

The Anti-STAT1 (Phospho-Ser727) and Anti-STAT1 Antibody were used to measure the effects of stimulation. Data in Figure 5 shows upregulation of the phosphorylated target protein detected via fluorometric reaction.

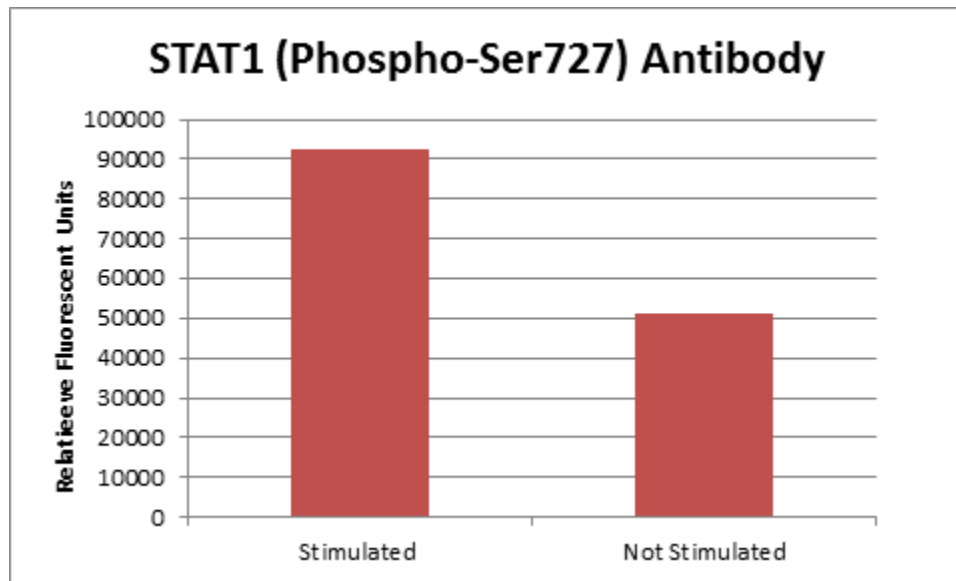

**\*\* Results may vary based on stimulations and fluorometer used\*\***

**Figure 5.** Hek 293 cells were starved for 24 hrs. and stimulated with PMA 200nM for 2 hrs.

## REAGENT PREPARATION

Ensure that all supplied solutions are at ambient temperature before use. It is recommended to conduct assays for all controls and samples in duplicate. Sufficient reagents are provided to assay 192 wells; therefore, on the day of the experiment, prepare only the required amount. During incubation steps, use an orbital shaker set at 200 rpm to ensure proper equilibration of solutions in the cell culture plate wells. **All other provided solutions are ready to use.**

1. Dilute the 10x TBS to 1x TBS by combining 9 volumes of ddH<sub>2</sub>O and 1 volume of 10x TBS.
2. Dilute the 15x Wash Buffer to 1x Wash Buffer using 14 volumes of ddH<sub>2</sub>O and 1 volume of 15x Wash Buffer.
3. Prepare the Fixing solution to 4% paraformaldehyde (for adherent cell lines) or 8% paraformaldehyde (for suspension or loosely attaching cell lines). Dilute 37% paraformaldehyde with 1x PBS to desired paraformaldehyde percentage.
4. For suspension or loosely attached cells: Prepare Poly L Lysine to working concentration as stated in manual. (Not included, recommended Sigma Cat# P4832).

*If the cell line is known to digest Poly-L-Lysine, Poly-D-Lysine may be an alternative for this step. Follow the manufacturer's instructions if using a different Poly-L-Lysine or Poly-D-Lysine.*

5. Immediately before use, add 50 ul of 100x Anti-STAT1 (Phospho-Ser727) Antibody and 50 ul of 100x Anti-GAPDH Antibody to 4,900 ul of Primary Antibody Diluent (for one plate, 96 wells). Gently mix and label the tube as "Primary Antibody Mixture Phospho".
6. Immediately before use, add 50 ul of 100x Anti-STAT1 Antibody and 50 ul of 100x Anti-GAPDH Antibody to 4,900 ul of Primary Antibody Diluent (for one plate, 96 wells). Gently mix and label the tube as "Primary Antibody Mixture Non-Phospho".
7. Immediately before use, for one plate, mix 3 ml of Dye-1 Conjugated Anti-Rabbit IgG Antibody and 3 ml of Dye-2 Conjugated Anti-Mouse IgG Antibody. Gently mix and label tube as "Secondary Antibody Mixture". This solution is light sensitive. Please store and handle in the dark.

## ASSAY PROTOCOL

1. Seed 100 µl of 20,000 adherent cells into each well in a 96-well plate.

**Note:** Optimal seeding concentrations should be determined by investigator and experimental goals.

For suspension cells and loosely attached cells:

- a. Coat plates evenly with 20 µL of Poly-L-Lysine (Sigma Cat# P4832).
  - b. After 5 minutes, remove solution by aspiration and rinse once with tissue culture grade water.
  - c. Allow 2 hours to dry prior to seeding cells.
2. Incubate cells overnight or at least 6 hours at 37°C, 5% CO<sub>2</sub>.
  3. Treat the cells as desired.
  4. Gently aspirate the cell culture medium and rinse each well twice with 200 µl of 1x TBS.

**Note:** For loosely attached cells, omit this 1x TBS rinsing step to avoid cell loss. Instead, proceed directly to the fixing instructions in Step 5 after removing the culture medium.

5. Gently remove the 1x TBS and fix the cells with 100 µl of 8% paraformaldehyde per well. Incubate for 20 minutes at room temperature (RT).

**Caution:** Paraformaldehyde is volatile. Use appropriate personal protective equipment (PPE) such as gloves, mask, and goggles.

6. Aspirate the solution and add 300 µl of 1x Wash Buffer to each well. Gently shake the plate on an orbital shaker for 2-3 minutes. Repeat this process 3 times. After the final wash, invert the plate and tap it against clean paper towels to remove any remaining liquid.

**Note:** After this step, plates can be stored at 4°C for up to a week.

7. Add 100 µl Quenching Buffer to each well and incubate for 20 mins at RT.
8. Repeat the washing process as described in Step 6.
9. Add 200 µl of Blocking Buffer and incubate for 1 hour at RT.
10. Repeat the washing process as described in Step 6

11. Add 50  $\mu$ l of 1x Primary Antibody Mixture solutions into each well and incubate overnight at 4°C. If target concentration is known to be high, incubate for 2 hours at RT. Designate 2 wells for negative control by omitting the 1x Primary Antibodies.
12. Repeat the washing process as described in Step 6.
13. Add 50  $\mu$ l of Secondary Antibody Mixture solutions into each well and incubate for 1.5 hrs.
14. Repeat the washing process as described in Step 6
15. Read the plate(s) at Ex/Em: 651/667 (Dye 1) and 495/521 (Dye 2). Shield plates from direct light exposure.

## DATA ANALYSIS

### Background subtraction

Average the duplicate or triplicate absorbance readings for each control and sample, subtracting them from the averaged absorbance for wells where primary antibodies have been omitted.

### Anti-GAPDH Internal positive control

The OD450 values obtained ensure the assay is functioning, and cell densities between wells should be proportionate to seeding concentration. Stimulation of cells should impact target primary antibody absorbance values, with no effect on the anti-GAPDH absorbance.

### Intra well Anti-GAPDH Normalization

The Anti-GAPDH normalization method enables intensity normalization within the same well. By using the ratio, the fluorescent Em.667 readings obtained can be normalized with the Em.521 values. This constitutes a within-well method of analysis.

a.  $\text{Em.667 anti-STAT1 (Phospho-Ser727)/Em.521 Anti-GAPDH.}$

b.  $\text{Em.667 anti-STAT1/Em.521 Anti-GAPDH.}$

By doing so, the expression levels of STAT1 (Phospho-Ser727) and STAT1 are adjusted to account for the cell density present in each well.

### Phosphorylation to Non-Phosphorylation Comparison

After Anti-GAPDH Normalization, the resultant proportional values can be used to analyze the effects of stimulants in this assay.

A non-phospho specific primary antibody is included to normalize the absorbance values between phosphorylated to non-phosphorylated values.

Both antibodies utilize the proportion for analysis as follows:

$\text{Em.667 (Anti-GAPDH Normalized) (STAT1 (Phospho-Ser727))}/$

$\text{Em.667 (Anti-GAPDH Normalized) (STAT1)}$

## **SHORT PROTOCOL**

**Seed cells into wells and incubate overnight at 37°C, 5% CO<sub>2</sub>**



**Apply desired treatment conditions**



**Add 100ul of Fixing Solution per well and incubate 20 minutes at room temperature**



**Add 100ul per well of Quenching Buffer and incubate 20 minutes at room temperature**



**Add 200ul per well of Blocking Buffer and incubate for 1 hour at room temperature**



**Add 50ul of “Primary Antibody Mixture P” and/or “Primary Antibody Mixture NP” per well and incubate overnight at 4°C**



**Add 50ul of “Secondary Antibody Mixture” per well and incubate for 1.5 hours at room temperature**



**Read the plate at Ex/Em: 651/667 (Dye 1) and 495/521 (Dye 2)**

## HEALTH AND SAFETY PRECAUTIONS

- Reagents provided in this kit may be harmful if ingested, inhaled or absorbed through the skin.
- Fixing Solution contains formaldehyde. Formaldehyde is known to be a highly toxic reagent. Personal protection is strongly recommended while working with this chemical.

## ASSAY RESTRICTIONS

- This ELISA kit is intended for research purposes only, NOT diagnostic or clinical procedures of any kind.
- Materials included in this kit should NOT be used past the expiration date on the kit label.
- Reagents or substrates included in this kit should NOT be mixed or substituted with reagents or substrates from any other kits.
- Variations in pipetting technique, washing technique, operator laboratory technique, kit age, incubation time or temperature may cause differences in binding affinity of the materials provided.
- The assay is designed to eliminate interference and background by other cellular macromolecules or factors present within any biological samples. However, the possibility of background noise cannot be fully excluded until all factors have been tested using the assay kit.
- Individual results may vary due to differences in technique, plasticware and water sources.

## TROUBLESHOOTING GUIDE

| Problem                | Cause                          | Solution                                                                                                                          |
|------------------------|--------------------------------|-----------------------------------------------------------------------------------------------------------------------------------|
| <b>Low Signal</b>      | Storage and Expiration         | Read Manual for storage condition and expiration.                                                                                 |
|                        | Antibody Dilution              | Dilute antibody 1:100 with Antibody Diluent.                                                                                      |
|                        | Low Cell Number                | Seed more cells.                                                                                                                  |
|                        | Cells Detach from Well-Bottoms | Add solutions slowly from the side wall of the wells. Use Poly-L-Lysine to coat the plate period seeding the cells, if necessary. |
| <b>High Background</b> | Washing                        | Remove Wash Buffer completely.                                                                                                    |
|                        | High Cell Number               | Reduce the amount of seeded cells.                                                                                                |
| <b>Variation</b>       | Pipetting                      | Check and/or calibrate pipettes.                                                                                                  |
|                        | Washing                        | Remove Wash Buffer completely.                                                                                                    |
|                        | Cells Detach from Well-Bottoms | Add solutions slowly from the side wall of the wells. Use Poly-L-Lysine to coat the plate period seeding the cells, if necessary. |

**ELISA PLATE TEMPLATE**

|   | 1 | 2 | 3 | 4 | 5 | 6 | 7 | 8 | 9 | 10 | 11 | 12 |
|---|---|---|---|---|---|---|---|---|---|----|----|----|
| A |   |   |   |   |   |   |   |   |   |    |    |    |
| B |   |   |   |   |   |   |   |   |   |    |    |    |
| C |   |   |   |   |   |   |   |   |   |    |    |    |
| D |   |   |   |   |   |   |   |   |   |    |    |    |
| E |   |   |   |   |   |   |   |   |   |    |    |    |
| F |   |   |   |   |   |   |   |   |   |    |    |    |
| G |   |   |   |   |   |   |   |   |   |    |    |    |
| H |   |   |   |   |   |   |   |   |   |    |    |    |

## TECHNICAL SUPPORT

For troubleshooting, information or assistance, please go online to [www.assaybiotech.com](http://www.assaybiotech.com) or contact us at:

Assay Biotechnology Company, Inc.  
2200 Ringwood Ave.  
San Jose, CA 95131  
United States of America

Email: [tech@assaybiotech.com](mailto:tech@assaybiotech.com)

Phone: (408) 747-0185  
Toll-Free Phone: (877) 883-7988

Fax: (408) 747-0145  
Toll-Free Fax: (877) 610-9758



**Over 3,000 ELISA Kits including Sandwich, Cell-Based and  
Transcription Factor ELISA Kits**

**Visit us at [www.assaybiotech.com](http://www.assaybiotech.com)**